

# Technical Specifications



- Overall dimensions \_\_\_\_\_ 2050 mm x 940 mm
- Mattress platform \_\_\_\_\_ 1800 mm x 720 mm
- Hi - Low \_\_\_\_\_ 900 mm - 580 mm
- Side rail height \_\_\_\_\_ 410 mm
- Side rail length \_\_\_\_\_ 1458 mm
- Backrest angle \_\_\_\_\_ 0 - 85°
- Safe working load \_\_\_\_\_ 220 kg
- Casters (Central locking) \_\_\_\_\_ 200 mm

**Technical and Quality Assurance Standards**  
 BS 5402  
 ISO 13485  
 ISO 14971

# Positions



**Piyatech Sdn. Bhd.**  
 B3-U5-7, Solaris Dutamas, Jalan Dutamas 1,  
 50480 Kuala Lumpur, Malaysia  
 Tel: (+60 3) 62113120, 62112120, 62112120  
 www.piyatech.com  
 email: info@piyatech.com




To improve products, Piyatech Sdn. Bhd. reserves the right to make changes without notice in design, specifications, and models.

# Piyatech

## Patient Trolley 8000T1 Type 1



Piyatech 2023-8000T1 type 1 V01 Validity Dec. 2024

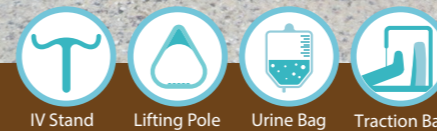

MDA Establishment License No.: MDA-1216-K120  
 MDA Product Reg. No.: GMD45614744918A



## Features



### Side Rails

8000T1 has a pair of durable, robust, and easy to clean side rails that provide a protective height and length to ensure the patient's safety during transportation and treatment period.

Both safety side rails are equipped with innovative and easy to use locking mechanisms that allow the rails to be securely locked in an upright position to maximize the patient's safety.

Pulling the upward locking handle, the caregiver can easily collapse the side rails.



### Casters

200 mm central locking casters run extremely smoothly, the quality that makes it significantly easy for the nursing staff to move the bed. In addition, the foot pedals allow effortless transition between the three modes of the central locking system:

- Safe central locking (down)
- Easy maneuvering and positioning (middle)
- Fast directional movement (up)

to provide a smooth and controlled single person transportation of the bed in the medical centers.



### Hydraulic Cylinder

The lower section of the 8000T1 bed is equipped with a foot pedals on each side, operated by robust and high-quality hydraulic cylinder. This pedal allow for easy adjustment of the trolley height.



### Adjustable Grip Handles

Two pairs of adjustable handles at the head and foot end of 8000T1 allow the nursing staff to move the trolley easily in the medical centers.



### Mattress Platforms

The trolley's ABS mattress platforms boast many helpful features, like smooth surfaces, high durability, and a lightweight, easily transportable frame enabling quick and painless removal for cleaning and disinfection, including steam cleaning.



### Oxygen Cylinder Holder

Oxygen cylinder holder is located underneath the trolley.



### Urine Bag Holder

Urine bag holders can easily install on either side of the trolley in one step.