

Technical Specifications



- Overall dimensions ————— 2100 mm × 920 mm
- Mattress platform ————— 1910 mm × 660 mm
- Hi-Low ————— 900 mm - 600 mm
- Side rail height ————— 310 mm
- Side rail length ————— 1620 mm
- Backrest angle ————— 0 - 39°
- Bed weight ————— 158 kg
- Safe working load ————— 220 kg
- Casters ————— 200 mm
- Trend. / Reverse Trend. ————— ± 10°
- Fifth Caster (Optional)

Technical and Quality Assurance Standards
 BS 5402
 ISO 13485
 ISO 14971

Positions



Piyatech Sdn. Bhd.
 B3-U5-7, Solaris Dutamas, Jalan Dutamas 1,
 50480 Kuala Lumpur, Malaysia
 Tel: (+60 3) 62113120, 62112120, 62112120
 www.piyatech.com
 email: info@piyatech.com




To improve products, Piyatech Sdn. Bhd. reserves the right to make changes without notice in design, specifications, and models.

Piyatech

Patient Trolley 8000T4



Piyatech 2023-8000T4 V01 Validity Dec. 2024

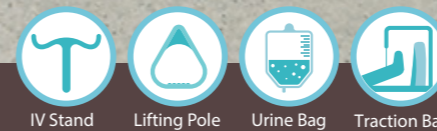

MDA Establishment License No.: MDA-1216-K120
 MDA Product Reg. No.: GMD45614744918A




Features




■ Side Rails

8000T4 has a pair of durable, robust, and easy to clean folding aluminum safety side rails that provide a protective height and length to ensure the patient's safety during transportation and treatment period.

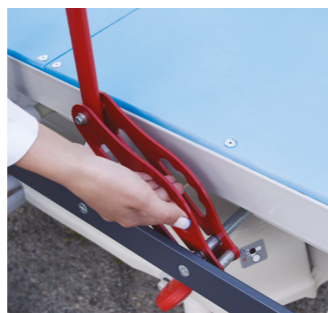

■ Gas Spring

With the aid of a pair of gas springs situated beneath the mattress platform in the head section, the caregiver can effortlessly and swiftly adjust the backrest angle manually.

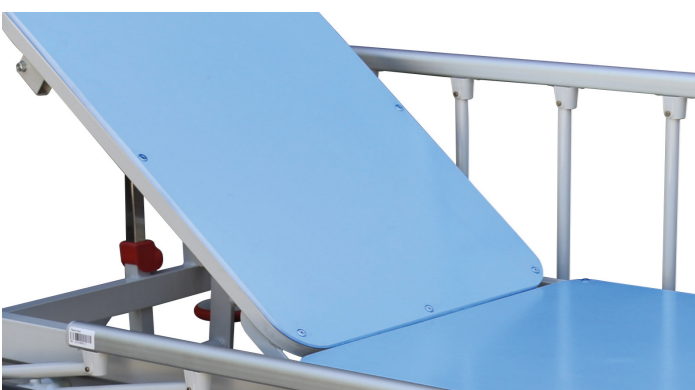


■ Mattress Platforms

The patient trolley's mattress platforms are made of durable phenolic board (compact board) with smooth surfaces, allowing for effortless cleaning and disinfection, while also being compatible with steam cleaning for enhanced hygiene.



■ Adjustable Grip Handles

Two pairs of adjustable handles located at either end of the 8000T4 allow the nursing staff to securely handle trolley during the transportation.

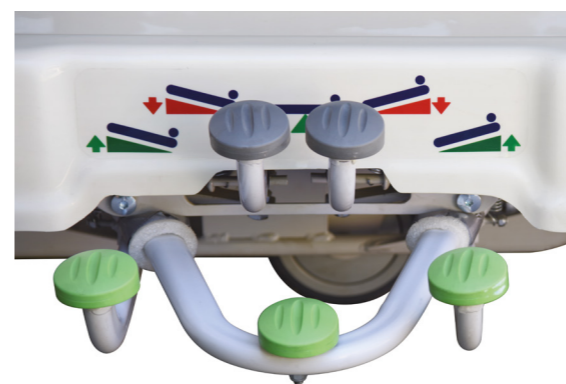

■ Hydraulic Cylinder

The lower section of the 8000T4 bed is equipped with five foot pedals on each side, operated by robust and high-quality hydraulic cylinders. These pedals allow for easy adjustment of the back, knee, high, low, Trendelenburg, and reverse positions.



■ Casters


200 mm central locking casters run extremely smoothly, the quality that makes it significantly easy for the nursing staff to move the bed. In addition, the foot pedals allow effortless transition between the three modes of the central locking system to provide a smooth and controlled single person transportation of the bed in the medical centers.

- safe central locking (down)
- easy maneuvering and positioning (middle)
- fast directional movement (up)



■ Bumpers

Innovative rotary bumpers which are located on four corners of beds can prevent to a large extent the damages that may cause collision of the bed to the surroundings during transportation.

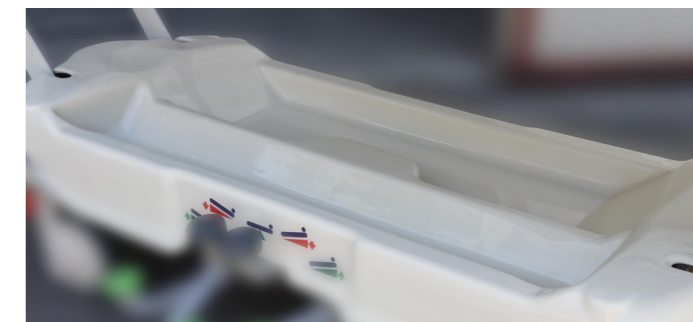

■ Oxygen Tank Holder and Utility Tray

The lower part of the trolley features a dedicated space for safe storage of the oxygen cylinder, accompanied by a utility tray for additional convenience.