

General Bed

25000MC1

Type 2




General Bed 25000MC1 Type 2

Technical Specifications



| | |
|--------------------|-------------------|
| Overall dimensions | 2260 mm x 1010 mm |
| Mattress platform | 2000 mm x 810 mm |
| Height | 440 mm |
| Side rail height | 410 mm |
| Side rail length | 1620 mm |
| Back rest angle | 0 - 70° |
| Knee rest angle | 0 - 35° |
| Bed weight | 90 kg |
| Safe working load | 220 kg |
| Casters | 125 mm |

Technical and Quality Assurance Standards
 IEC 60601-2-52
 ISO 13485
 ISO 14971

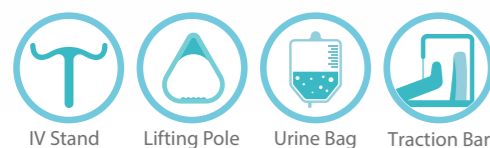

Features



Side Rails
 The parallel side rails of the 25000MC1 offer a unique and stylish design produced from a combination of blow molded plastic and steel bars to provide a protective height and length. The side rails adhere to the highest standards of safety to prevent patient falls and entrapment.



Side Rails Locking Mechanism
 Pushing up the red lever on either side of the bed engages the release mechanism allowing the side rails to fold down.



Foldable Cranks
 Back and knee positions can be adjusted by two foldable cranks located at foot end of the bed.

General Bed 25000MC1Type 2



Bed Ends

The seamless blow-molded bed ends are lightweight, durable, and easy to clean. The bed ends are of high quality and compatible with steam cleaning. To have better access to the patient's head in emergency cases, the caregiver can easily remove the bed ends by grasping the hand holds and lifting them up in one single step. Moreover, the bed ends have integrated steel fittings which make them extremely strong so the nursing staff can move the bed by holding and pushing them without worries.



Mattress Platforms

The injection-molded mattress platforms exhibit exceptional qualities, including smooth surfaces, high durability, and effortless removal for easy cleaning and disinfection. Additionally, the mattress platforms are designed to be compatible with steam cleaning, ensuring thorough and effective sanitation.



Casters

125 mm anti-static casters run extremely smoothly; the quality that makes it significantly easy for the nursing staff to transport the bed in the medical centers. Having diagonal locks, the casters can be locked to prevent the bed from unintended movement. (150 mm casters are available as an option)



IV Pole Holder

Four IV Pole Holders are Located on the head and foot end.



Urine Bag Holder

Beds are equipped with a urine bag holder which can be easily installed on either side of the bed.



Bumpers

Innovative rotary bumpers which are located on four corners of beds can prevent to a large extent the damages that may cause collision of the bed to the surroundings during transportation.

General Bed 25000MC1 Type 1



Safety Side Rail



Bed Ends

The seamless blow-molded bed ends that are lightweight, durable, and effortless to clean. These high-quality bed ends are compatible with steam cleaning. For enhanced accessibility to the patient's head during emergency situations, caregivers can conveniently remove the bed ends in a single step by grasping the hand holds and lifting them up. Additionally, the integrated steel fittings within the bed ends provide exceptional rigidity, enabling nursing staff to move the bed with ease and confidence.



Bumpers

The bed is equipped with bumpers positioned on all four corners, effectively mitigating potential damages caused by inadvertent collisions with the surroundings during transportation.