

Patient Trolley

**8000T4**

Type 2




# Patient Trolley 8000T4 Type 2

## Technical Specifications



Overall dimensions	2100 mm x 920 mm
Mattress platform	1910 mm x 660 mm
Hi-Low	900 mm - 600 mm
Side rail height	310 mm
Side rail length	1620 mm
Back rest angle	0 - 85°
Knee rest angle	0 - 45°
Bed weight	158 kg
Safe working load	220 kg
Casters	200 mm
Trend. / Reverse Trend.	± 10°
IP rating	IP66
Fifth Caster (Optional)	

**Technical and Quality Assurance Standards**  
 IEC 60601-1  
 BS 5402  
 ISO 13485  
 ISO 14971

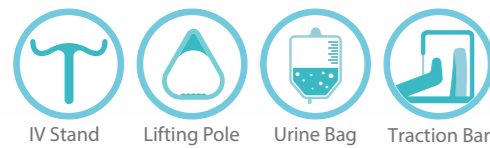

IV Stand   Lifting Pole   Urine Bag   Traction Bar



Back rest 0-85°   Knee rest 0-45°   Max. height 900 mm   Min. height 600 mm   Trendelenburg -10°   Reverse Trend. +10°

## Features

### Side Rails

8000T4 has a pair of durable, robust, and easy to clean folding aluminum safety side rails that provide a protective height and length to ensure the patient's safety during transportation and treatment period.



### Side Rails Locking Mechanisms

Both safety side rails are equipped with innovative and easy to use locking mechanisms that allow the rails to be securely locked in an upright position to maximize the patient's safety. Pulling upward locking handle, the caregiver can easily collapse the safety side rails.



### Handset Control

User-friendly handset controller with illustrative schematics makes it simple and easy for the caregiver to adjust the bed with different positions. The handset controller can be placed on both safety sides for easy access.



# Patient Trolley 8000T4



## Actuators

All positions of this bed are achieved by powerful, reliable, and low noise actuators. As the actuators have IP66 rating, they are protected against ingress of dust and water. This prominent feature of actuators makes it easy to clean and disinfect the bed in the medical centers on a regular basis.



## Mattress Platforms

The trolley's ABS mattress platforms boast many helpful features, like smooth surfaces, high durability, and a lightweight, easily transportable frame enabling quick and painless removal for cleaning and disinfection, including steam cleaning.



## CPR Mechanism

Two CPR quick release handles are located on either side of the trolley for easy and fast access in the case of an emergency. The manual quick-release mechanism operates independently from the electrical CPR of the trolley that makes it fully operational in spite of a power outage.



## Adjustable Grip Handles

Two pairs of adjustable handles located at either end of the 8000T4 allow the nursing staff to securely handle trolley during the transportation.



## Oxygen Tank Holder and Utility Tray

The lower part of the trolley features a dedicated space for safe storage of the oxygen cylinder, accompanied by a utility tray for additional convenience.



## IV Pole Holder

Four IV pole holders are located in each corner of the trolley.



## Casters

200 mm anti-static central locking casters run extremely smoothly, the quality that makes it significantly easy for the nursing staff to move the bed. In addition, the foot pedals allow effortless transition between the three modes of the central locking system to provide a smooth and controlled single person transportation of the bed in the medical centers.

- safe central locking (down)
- easy maneuvering and positioning (middle)
- fast directional movement (up)



## Bumpers

Innovative rotary bumpers which are located on four corners of beds can prevent to a large extent the damages that may cause collision of the bed to the surroundings during transportation.



## Urine Bag Holder

Beds are equipped with a urine bag holder which can be easily installed on either side of the bed.